

ProChimia Surfaces Sp. z o.o.

81-823 Sopot, ul. Zacisze 2

REGON: 220121035

KRS: 0000242574

NIP: 5851419173


www.prochimia.com

Sopot, 9.03.2015


Custom PDMS STAMP Symbol: MP-SQ 5x5 micron

Description: The MP-SQ 5x5 micron PDMS stamp is 1 cm square size with protruding square features of 5 microns (height of 4 microns) arranged in symmetrical array of rows and columns with 5 microns period between neighboring features.

General scheme of PDMS stamp pattern is presented on Fig.1

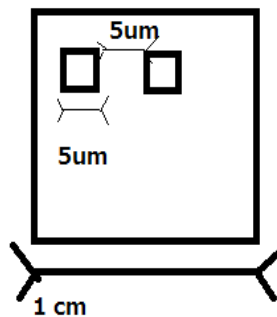

PDMS Stamp MP-SQ 5x5 micron image (Camera: Nikon LV100ND)

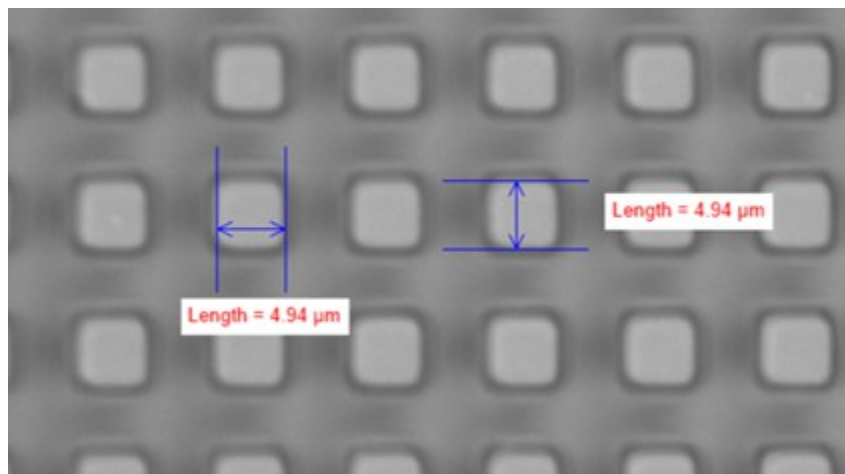

Storage: Store in room temperature of + (4-8°C) in dust-free chamber (original packaging or petri dish).

Cleaning: For cleaning rinse with absolute ethanol and dry in a stream of nitrogen. For any dirt removal submerge in ethanol and sonicate for 5 min. in ultra wave bath. Subsequently dry in a stream of nitrogen.
